
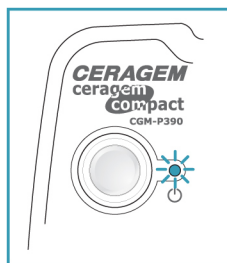



## 4.Using Your Product

- 1 Turn on the power switch on the side of the product.

After turning on the power switch, the Power LED will flash until the product starts operating.

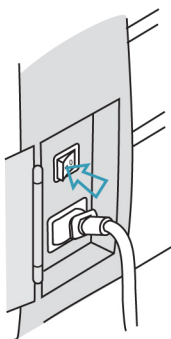


- 2 Press the Power Button on the front of the product. (Manual Mode)

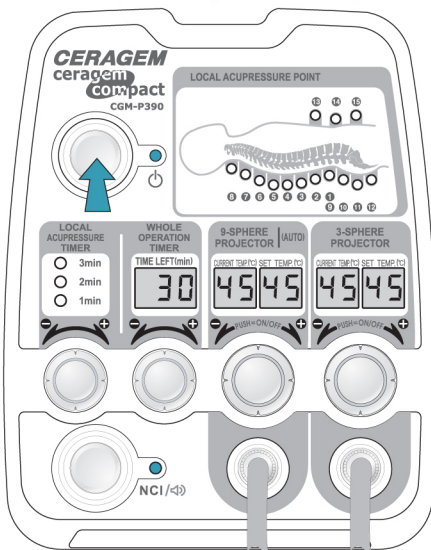

- 3 **Negatively Charged Ion (NCI) Function**

Press the NCI button to turn the NCI function ON or OFF.

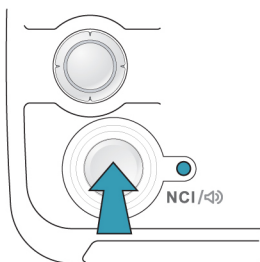

### Audible Sound On/Off

Press the Audible Sound button for 1~2 seconds to turn the Audible Sound function ON or OFF.

- 4 Set the temperature to the desired level by turning the temperature control for the corresponding Projectors. Temp. Range: 30-60 °C

### Turning Each Projector ON or OFF


Briefly press the temperature controller for the corresponding Projector.

## 5 Auto/Manual Mode

### Auto Mode

Press the 9-Sphere Projector Temperature Controller for 1~2 seconds.

(See Page 10 for operation instructions.)

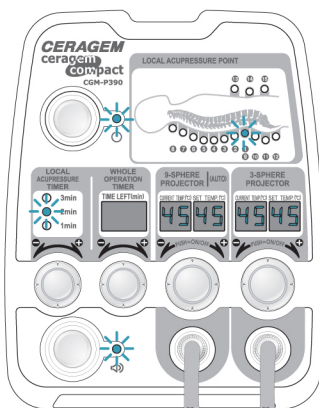


< An example of the auto mode. >

The local acupressure point display and the local acupressure timer are illuminated.

### Manual Mode

- Turn the power button on and the manual mode will operate.
- Press and hold the 9-Sphere projector button for 1~2 seconds in the auto mode to switch to the manual mode.
- Turn the 9-Sphere projector off in the auto mode will be transfer to the manual mode.
- Upon completion of the auto mode, the manual mode will operate.
- Set the Operation Time: 1 ~ 120 minutes
- The operation stops when the time display reaches [0].

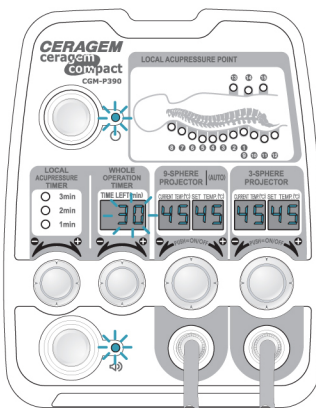

< An example of the manual mode. >

The whole operation timer is illuminated.

- 6** To operate again, press the power button. When you are finished, switch off the power switch on the side of the product.

## ● Messages

| Display                                                                                     | Cause & Solution                                                                                                                                                                                                                                                                   |
|---------------------------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| <p>-L</p>  | <p>① When the Projectors are not connected, the [-L] message is displayed on the FND window.</p> <p>② If the temperature is lower than -10°C "-L" is shown on the display.</p> <p>③ If the temperature is between -10°C and 0°C "0" is shown on the display.</p>                   |
| <p>E1</p>  | <p>When the Projectors are not connected properly, the [-L] message will be displayed on the FND window, then it will be switched to the flickered [E1] mode for 5 seconds, then the power automatically turns off.</p> <p>※ Connect the Projectors.</p>                           |
| <p>E2</p>  | <p>When a Projector has over-heated, the [E2] message is displayed on the FND window (Entire Operation Timer) for 5 seconds, then the power automatically turns off.</p> <p>※ Do not dismantle the Product.</p> <p>Contact the Ceragem Retailer or Service Center nearest you.</p> |